

Guidance Document

EN ISO 15223-1 new symbols for SBS

The Sterile Barrier Association (SBA) created and validated new symbols for Sterile Barrier System (SBS) configurations respective packaging systems for sterile medical devices for inclusion into ISO 15223-1. Reasons for inclusion of such symbols are to control specific risks with aseptic presentation, to comply with new legal requirements deriving from the EU-MDR 2017/745 and to provide additional user benefits.

Why do we need symbols for identification of SBS configurations?

- Sterile medical products are manufactured and processed by Medical Device Manufacturers (MDM's) and in healthcare or reprocessing facilities.
- Sterile packaging systems prevent ingress of microorganisms and allow aseptic presentation.
- Sterile packaging systems are composed of at least one sterile barrier system which maintains sterility and allows for aseptic presentation
- A Protective Packaging (PP) layer is often added to physically protect the SBS and its contents until the point of use. Protective packaging can be outside and also inside the SBS.
- In many cases, there is no difficulty differentiating the two. A corrugated outer dispenser box is obviously not a sterile barrier system but may still be designed to be used as protective packaging for transport and storage. A single sterile barrier system, e.g. a pouch, containing a sterile medical product, may be easily identified as a sterile barrier system.
- There are circumstances however, where it is difficult to differentiate between a validated sterile barrier system and protective packaging that looks like a sterile barrier system. In these cases, risks could arise during aseptic presentation: risks of contaminating the device and/or the sterile field and/or sterile gowns of operating room personnel

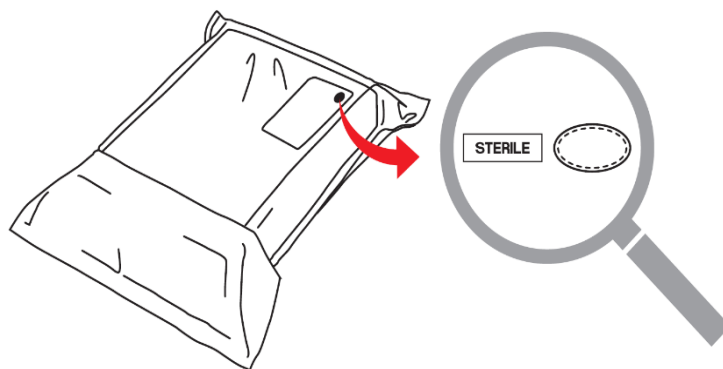
Symbols help to identify and differentiate between SBS and PP to mitigate that risk:

The new symbols are composed from ovals, which are formed either from:

- a solid line which indicates a Sterile Barrier System layer (maintaining sterility) or
- a dashed line which indicates a Protective Packaging layer that is not a validated microbial barrier.

The symbols shall be printed on the label which identifies the medical device, adjacent to or in combination with the symbol 'sterile'. A typical packaging system configuration for sterile medical products could be made from a header bag, containing sterile products which are wrapped in protective packaging. The protective packaging does not provide validated barrier properties but is used to provide an aseptic presentation tool. In the diagram shown below:

- The outer solid line oval represents the header bag which is the Sterile Barrier System.



- The dashed line oval indicates that the medical product is wrapped in a protective packaging layer for aseptic presentation.

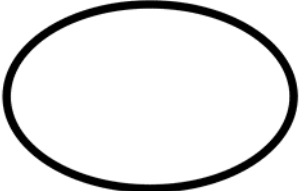
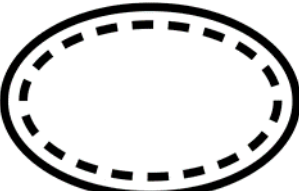
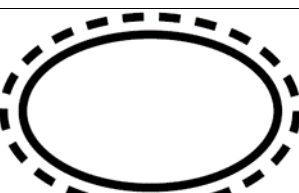
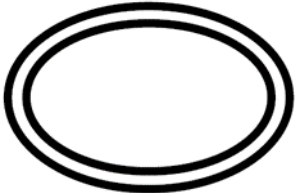

- The protective packaging layer does not provide validated barrier properties.

- Operating room staff can easily identify this particular packaging configuration, indicating they should place only the inner pack on the sterile field for aseptic presentation.

Guidance Document

EN ISO 15223-1 new symbols for SBS

All new symbols in an overview:

Symbol	What it represents	Recommended handling / usability
	Ref: ISO 7000 - 3707 Single sterile barrier system	Aseptic presentation technique requires opening by an assistant nurse. (Sterile) Scrub nurses or surgeons must not touch the outer surface of the packaging. Pack must not be placed on sterile surfaces.
	Ref: ISO 7000 - 3708 Single sterile barrier system with protective packaging inside	Aseptic presentation technique requires opening of the outer packaging by an assistant nurse. Sterile nurses or surgeons must not touch the surface of the outer packaging. The inner layer with the sterile product may be handled by sterile personnel. Product in inner layer can be placed on sterile surfaces.
	Ref: ISO 7000 - 3709 Single sterile barrier system with protective packaging outside	Aseptic presentation technique requires opening by an assistant nurse. Sterile nurses or surgeons must not touch the outer surface of the sterile packaging. Pack must not be placed on sterile surfaces.
	Ref: ISO 7000 - 3704 Double sterile barrier system	Aseptic presentation technique requires opening of the outer sterile packaging by an assistant nurse. Sterile nurses or surgeons must not touch the outer surface of the sterile packaging. Outer packaging must not be placed on sterile surfaces. The inner sterile packaging may be handled by sterile personnel and can be placed on sterile surfaces.
	Double sterile barrier system with protective packaging outside	Aseptic presentation technique requires opening of the outer protective packaging and the outer sterile packaging by an assistant nurse. Sterile nurses or surgeons must not touch the outer surface of the sterile packaging. Outer packaging must not be placed on sterile surfaces. The inner sterile packaging may be handled by sterile personnel and can be placed on sterile surfaces.

Note: The 3-layer symbol has been validated but will not be part of ISO 7000 nor of ISO 15223-1.

Detailed descriptions on interpretation and application of the symbols are available in [ISO 15223-1:2016](#) and [ISO DIS 15223-1:2020](#)

With the exception of the 3-layer version, all symbols are available as graphical data files from the ISO 7000 online shop.

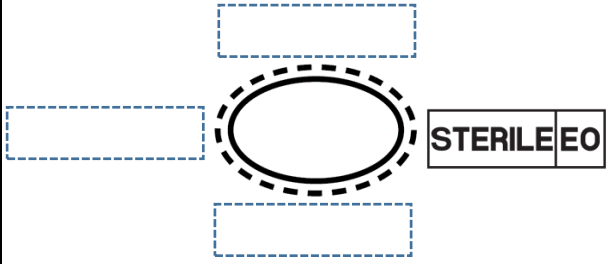

<https://www.iso.org/obp/ui#iso:grs:7000:3704> <https://www.iso.org/obp/ui#iso:grs:7000:3707>
<https://www.iso.org/obp/ui#iso:grs:7000:3708> <https://www.iso.org/obp/ui#iso:grs:7000:3709>

Guidance Document

EN ISO 15223-1 new symbols for SBS

The placement of the symbols on the packaging:

According to ISO DIS 15223-1:2020, the symbols shall be placed on the label which identifies the medical device, adjacent to or in combination with the symbol 'sterile'. The following combinations are possible:

Placement next to the symbol "sterile"	Placement within the symbol "sterile"
 <p>Depending on the overall label design</p>	 <p>In case of small packages</p>

Notes

- The symbol sterile will indicate the method of sterilisation, the table above indicates a medical device that has been sterilised using ethylene oxide as an example
- The manufacturer should determine the appropriate size for the symbol to be legible for its intended use.

Which layers of packaging should be labelled?

The MDR requires the packaging which maintains the sterile condition of a device ('sterile packaging') to have an indication permitting the sterile packaging to be recognised as such. The MDR does not specify which layers of packaging must be labelled.

ISO 11607-1: 2019 requires in subclause 6.1.8: "If the packaging system to be opened at the point of use consists of more than one packaging layer, the sterile barrier system(s) shall have an indication to be recognised as such."

According to ISO 11607, there is no requirement to label anything other than the SBS. However, the MDR also requires in annex 1 - GSPR 11.1, that the design allows for easy and safe handling and [...] prevent microbial contamination to eliminate, or reduce as far as possible, the risk of infection to the patient.

The decision to label protective packaging should be an outcome of the risk evaluation, and packaging system design process in order to achieve acceptable usability for aseptic presentation. If protective packaging looks like an SBS, then the validated symbols should be used to control the risk of unintentional contamination of the sterile field. In the case of paper board protective packaging, the risk is low that healthcare professionals consider this an SBS, so no labelling would be required. Depending on the user requirements and the intended use, the manufacturer can decide the need to label the protective packaging if it is considered useful for healthcare professionals.

The objective of the usability evaluation for aseptic presentation, a new requirement of ISO 11607-1: 2019, is to demonstrate evidence that the design including the respective labels allows for easy and safe handling.

Guidance Document

EN ISO 15223-1 new symbols for SBS

About the SBA

The Sterile Barrier Association (SBA) is the European trade association for companies who produce Sterile Barrier Systems (SBS) and associated equipment and accessories for the healthcare industry. Its mission is to be the recognised expert association in the healthcare industry, promoting the use of and providing education on the most suitable single use sterile barrier systems to ensure patient safety.

Most of the SBA members manufacture in Europe, many are global companies. All members are registered to ISO 9000 or another recognised higher level quality management system and many incorporate elements of GMP in their protocols. The majority are certified to EMAS or ISO 14001 as an environmental management system.

Director General

June 2020